

WELCOME

Congratulations on your new Hypervolt GO

You are one step closer to moving better, feeling better, and being a better version of you. The Hypervolt GO was created to keep you doing what you love. Whether you're an elite athlete, weekend warrior, or simply trying to get your body to feel its best, the Hypervolt GO will help you get there... and then some.

NEXT STEPS

Let's get started

1 Register your device at hyperice.com/register-product to activate your warranty and ensure easy returns, repairs, or refunds.

2 Download the **Hyperice App** to get the most out of your new Hypervolt GO, or any Hyperice device.



3 Charge your device by plugging in your new Hypervolt GO with the provided Hyperice charger. Fully charge the battery for up to six hours before first use.

Have questions? Call 1.855.734.7224, chat live at hyperice.com, or email: customersupport@hyperice.com.

Moving better, made simple

1 Power up

With your hand on the handle, press down on the power/speed button with your index finger. Hold for two seconds until the LED band (bottom of handle) illuminates and release. **Your Hypervolt GO is now ready for use.**

2 Adjust speeds / power off

Choose your desired speed (1-3) by pressing the power/speed button once per speed. The blue lights on the bottom of the handle will indicate speed setting. To power off, simply hold the power/speed button down for two seconds.

QUICK START GUIDE

3 Changing attachments

Choose desired head attachment and insert directly into the shaft/opening while pushing firmly and slightly twisting.

4 Recharging


Only use the provided Hyperice wall charger to uphold warranty. Full charge is indicated when the LED band on the handle stops flashing and becomes solid green.

5 Caring for your 'GO'


Make sure the power is OFF and the battery charger is not attached. Use a damp, clean, cloth and gently wipe your Hypervolt GO. Store in a clean, cool place when not in use.

HEAD ATTACHMENTS

Choosing the best head attachment

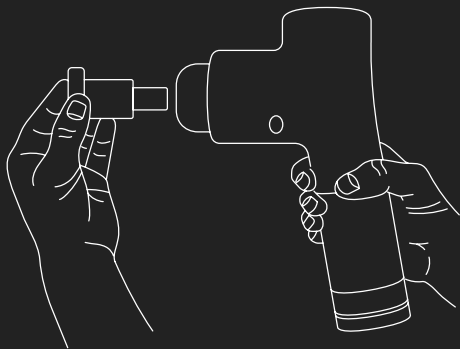

HEAD ATTACHMENTS

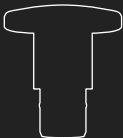

Flat

Larger surface areas such as quads, hamstrings, back, and chest.



Bullet


Trigger points found in glutes, hips, traps, shoulders, and other complex muscle groups.

Want to learn more? Visit hyperice.com/how-to for video tutorials on your Hypervolt GO or any Hyperice device.

HYPERICE APP

The next generation of Hyperice

Transform the way you move with the latest innovation by Hyperice. Powered by HyperSmart™, the Hyperice App seamlessly blends your physical and digital activity to create an experience entirely customized to you.

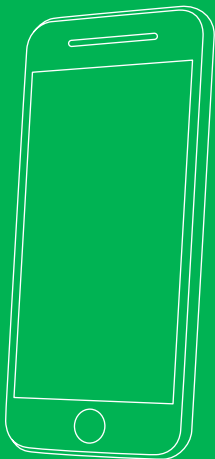

HYPERICE APP



Connect to health and wellness apps

Connect to apps like Strava and Apple Health for a more personalized experience.



Guided sessions

Guided warmup, recovery, and maintenance sessions.



Exclusive content

Exclusive content from the world's best athletes, trainers, and physiotherapists.



Operate your devices

Operate your Hyperice Bluetooth® connected devices.

Apple Health is a registered trademark of Apple, Inc. and used under license. Strava is a registered trademark of Strava, Inc. and is used under license. The Bluetooth® word mark and logos are a registered trademark owned by Bluetooth SIG, Inc. and is used under license.

HYPERCARE

A dedicated team just for you



Every Hyperice product is backed by a 1-year warranty and award-winning support from our Hypercare™ team – a group of specialists dedicated to your overall well-being and expert guidance on Hyperice products.

Call 1.855.734.7224

Chat [hyperice.com](https://www.hyperice.com)

Email customersupport@hyperice.com

HYPERCARE



Product knowledge

Answer questions about your device and best practices.



Place / track new orders

Help with placing new orders and working with shipping carriers.



Troubleshoot

Help troubleshoot and solve issues with your device.



Billing solutions

Research and provide solutions for any billing issues.




Initiate returns

Provide instructions for returns, refunds, or repairs.


OTHER PRODUCTS

A full suite of products to help you move better

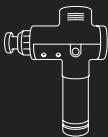

Hypervolt / Hypervolt Plus

Handheld Percussion
Massage Device

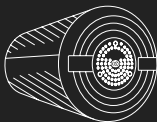


Vyper 2.0

High-Intensity
Vibrating Fitness
Roller




Hypersphere

Vibrating Massage Ball

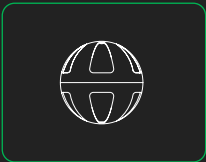


Hypersphere Mini

Compact Vibrating
Massage Ball

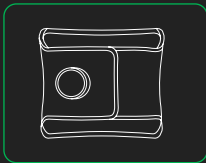

Venom

Wearable Heat +
Vibration

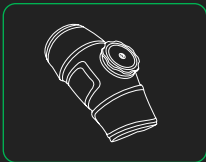

ICT

Ice Compression
Technology

Learn more / shop @ hyperice.com