Hypervolt Go 3

Hypervolt Go 3 筋膜枪
Hypervolt Go 3 筋膜槍
Hypervolt Go 3 Masajeador
Model: HVG3001

EN Operating Instructions

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EN

TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, OR PROPERTY DAMAGE, THIS DEVICE MUST BE USED IN ACCORDANCE WITH THE FOLLOWING WARNINGS, CAUTIONS, AND SAFETY INSTRUCTIONS.

IMPORTANT SAFETY INSTRUCTIONS - ORIGINAL INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following. Read all instructions before using the Hypervolt Go 3.

DANGER

To reduce the risk of electric shock:

- Always unplug this device from the electrical outlet immediately after using and before cleaning.
- Do not reach for a device that has fallen into water. Unplug immediately.
- Do not use while bathing or in a shower.
- Do not place or store device where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

WARNING

To reduce the risk of burns, fire, electric shock, or injury to persons:

- A device should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- This device can be used by children aged from 8 years and above and persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Use this device only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this device if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if it has gotten wet. Return the device to a service center for examination and repair.
- Do not carry this device by supply cord or use cord as a handle.
- Keep the cord away from heated surfaces.
- Never operate the device with the air openings blocked. Keep the air openings free of lint, hair, or any substance which may impair the flow of air.
- Never drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) devices are being used or where oxygen is being administered.
- Never operate on a soft surface such as a bed or couch where the air openings

may be blocked.

- Do not use device in close proximity to loose clothing or jewelry.
- Keep long hair away from device while in use.
- Cleaning and user maintenance shall not be done by children without supervision.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. After one hour use, allow device to rest for 30 minutes before re-using.
- Unplug the device after charging or prior to use.
- Use device only according to Operating Instructions.
- When using the Hypervolt Go 3 use only on dry, clean surfaces of the body by lightly pressing and moving across the body for approximately sixty(60) seconds per region. The Hypervolt Go 3 should only be used on soft tissue areas of the body and should not be used on any hard or bony areas of the body, including the head. Please immediately stop use of the Hypervolt Go 3 upon experiencing any pain or discomfort, apart from light muscle soreness. Bruising can occur from use of the Hypervolt Go 3, regardless of the pressure setting, and if bruising does occur you should not use the Hypervolt Go 3 on the bruised area until the bruise is fully healed. Use of the Hypervolt Go 3 should also be avoided on all bruises, contusions, rashes or irritated or injured areas of the skin until they are fully healed. When using the Hypervolt Go 3 be sure to keep fingers, toes, hair, and other body parts away from the back of the attachment head to avoid pinching or hair entanglement.
- The Hypervolt Go 3 should only be used on muscle tissue areas. The Hypervolt Go 3 should not be used on any hard or bony areas of the body, including the head and spine, the sides or front of the neck, the inner groin, or the back of the knee. Please immediately stop use of the Hypervolt Go 3 upon experiencing any pain or discomfort, apart from light muscle soreness. Bruising can occur from use of the Hypervolt Go 3 regardless of the pressure setting, and if bruising does occur you should not use the Hypervolt Go 3 on the bruised area until the bruise is fully healed. Use of the Hypervolt Go 3 should also be avoided on all bruises, contusions, rashes or irritated or injured areas of the skin until they are fully healed. When using the Hypervolt Go 3 be sure to keep fingers, toes, hair, and other body parts away from the back of the attachment head to avoid pinching or hair entanglement.

The Hypervolt Go 3 is intended to be charged with the Hyperice USB-C cord provided or by a 5V USB. Do not charge the Hypervolt Go 3 overnight or leave the Hypervolt Go 3 unattended when charging or in use. Use only with Listed/Certified ITE Power Supply. Do not use unapproved power source.

Please do not use the Hypervolt Go 3, or any percussion device without first obtaining approval from your physician if any of the following apply:

- Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pacemakers, recent surgery or injury, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any percussion or vibration device. These contraindications do not mean that you are not able to use a percussion or vibration device but we advise that before doing so you consult your physician. Ongoing research is being done on the effects of percussion massage for specific medical disorders that may result in the shortening of the list of contraindications as shown above.
- Use appliances only with specifically designated battery packs. Use of any other battery packs may create a risk of injury and fire.
- Under abusive conditions, liquid may be ejected from the battery; avoid contact. If contact accidentally occurs, flush with water. If liquid contacts eyes, additionally seek medical help. Liquid ejected from the battery may cause irritation or burns.
- Do not use a battery pack or appliance that is damaged or modified. Damaged or modified batteries may exhibit unpredictable behavior resulting in fire, explosion or risk of injury.
- Do not expose an appliance to fire or excessive temperature. Exposure to fire or

temperature above 130°C may cause explosion.

- Follow all charging instructions and do not charge the appliance outside of the temperature range specified in the instructions. Charging improperly or at temperatures outside of the specified range may damage the battery and increase the risk of fire.
- Have servicing performed by a qualified repair person using only identical replacement parts. This will ensure that the safety of the product is maintained.
- Do not modify or attempt to repair the appliance except as indicated in the instructions for use and care.

The Hypervolt Go 3 contains a battery which must be complied with safely at an appropriate e-waste disposal or recycling facility.

CAUTION

For household and commercial use. The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not disassemble, heat above 100° C or incinerate. Dispose of used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire.

SAVE THESE INSTRUCTIONS

The Hypervolt Go 3 is a handheld percussion massage device that delivers targeted pulses of pressure to care for muscles, relieve tension, provide a relaxing massage, accelerate warm up and recovery, and help maintain flexibility and range of motion.

Operating

Turn your Hypervolt Go 3 ON by pressing and holding the power setting button for 0.5 seconds. Once your device is powered ON the battery level indicator (LED band) at the base of the handle will illuminate. Turn the encoder dial to change speeds (1-5). The lights will indicate the current speed. To power off, simply hold the power button down for 0.5 seconds.

Changing head attachments

Only change attachments after speed setting is turned to off and head attachment stopped moving. To change head attachments, pull the existing head attachment straight out to remove and insert the desired head directly into the opening while pushing firmly.

Charging

Fully charge battery for up to four hours before first use. To charge, connect the USB-C cord into the charging port on the bottom of the handle and plug the cord into a USB-A adapter. Then plug the adapter into the wall. The LED light band will pulsate to show battery level and indicate active charging. The LED light band colors correspond to the charge level ranging from red (low) to green (fully charged). Full charge is indicated when the LED band remains green and fully illuminated. Battery may be recharged at any time and at any battery level. It is not recommended to fully deplete the battery to the red LED level. The average runtime is 2+ hours depending on speed level and pressure applied during use. It is recommended to turn the device OFF while the device is not in use or during storage and travel.


Cleaning and storage

Make sure the power is OFF and the battery charger is not attached. Use a damp,

clean cloth and gently wipe your Hypervolt Go 3. Store in a clean, cool place, out of direct sunlight, when not in use.

Specifications

1. Interchangeable head attachments (2)
2. Digital speed dial
3. Speed indicators (5)
4. ON / OFF power button
5. Battery level indicator (LED band)
6. Charging port

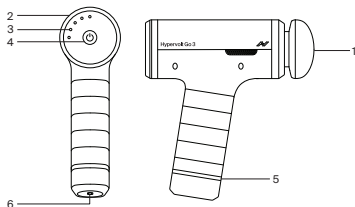5V --- 2A

Operating temperature: +32° F to +113° F [0° C to +40° C]

Storage temperature: -4° F to +158° F [-20° C to +60° C]

Charging temperature: +5°C to +40°C [+41°F to +104°F]

Rated Capacity: 7.2V 4250 mAh 30.6Wh

Frequency: Level 1 - 1500/25 Hz, Level 2 - 1750/29 Hz, Level 3 - 2000/33 Hz, Level 4 - 2250/38Hz, Level 5 - 2500/42Hz

Weight: 1.6 lbs/0.7 kg

Warranty

This product is covered by a limited warranty from Hyperice. Please visit hyperice.com/warranty to review the warranty in your country

Recycling

Please visit <https://hyperice.com/recycling> for information about recycling this packaging.

FOR MEXICO ONLY

- This device is not intended for use by persons (including children) whose physical, sensory or mental capabilities are different or reduced, or who lack experience or knowledge, unless such persons have been given supervision or training in the operation of the device by a person responsible for their safety.
- Children should be supervised to ensure that they do not use the devices as toys.
- The device must only be powered at the safety extra-low voltage indicated on the device.
- The battery must be removed from the device before disposal. The device must be disconnected from the charger when the battery is removed. The battery is removed for safe disposal.
- Operation of this equipment is subject to the following two conditions: (1) this equipment or device may not cause harmful interference, and (2) this equipment or device must accept any interference, including interference that may cause undesired operation.

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hyperice.com

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Hyperice, Inc., 525 Technology Drive, Suite 100, Irvine, California 92618 USA
Copyright 2025 © Hyperice, Inc. All rights reserved.
460-00131 Rev 3